




Diagnosis

Explorer	2
Excavator	5
Dressing Tweezers	6
Mouth Mirror	8
Disposable Instrument	10
Temporary Filling Instrument	10
Evacuation Tip	11
Stainless Dish / Tray	12



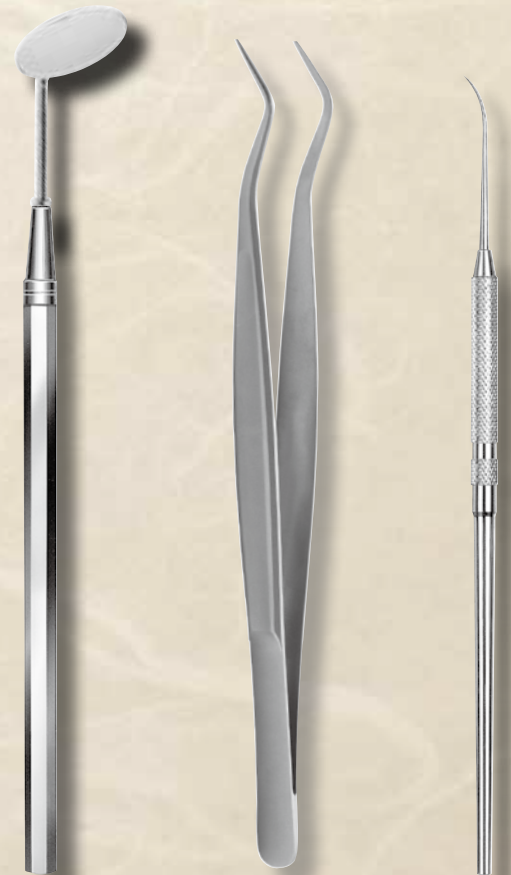

Explorer

-Single End-



#3
No. 01003



#5
No. 01005




#9
No. 01009



#23
No. 01023



#25
No. 01025

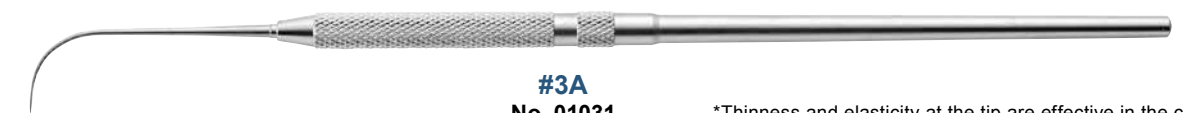

#3A
No. 01031

*Thinness and elasticity at the tip are effective in the check of dental tartar, cavity and surplus cement.



Multi
No. 01030


*A tip part is made thick for intensity.



#23BS
No. 01024



Octagonal Handle #9
No. 01109

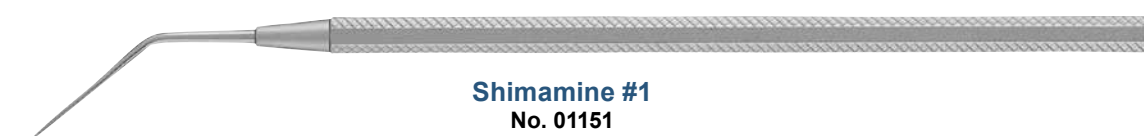


Shimamine #1
No. 01151

-Ball Point-




#1
No. 01221

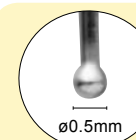

*The round tip realizes checking of caries without damaging pit and fissure tissues.
*Suitable for filling of sealant.


-Double End-

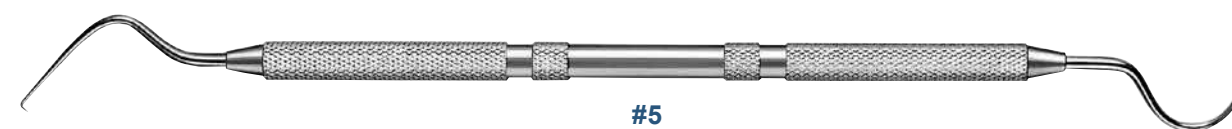

#5
No. 01205



#6
No. 01206

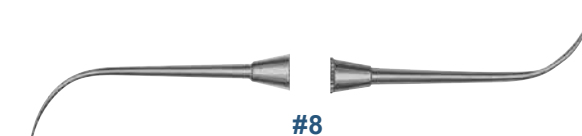

#8
No. 01208



#3-9
No. 01209



#11-12
No. 09224



*Explorer with the shank of "Gracey Curette #11-12".
*Thinness and elasticity at the tip are effective in the check of dental tartar, cavity and surplus cement.

- for Root Canal Therapy -



#16
No. 01316





NC #18
No. 11305

Size / 160mm



NC #18 with Serration
No. 11306

Size / 160mm



Mini #18
No. 11308

Size / 145mm




Mini #18 with Serration
No. 11309

Size / 145mm



NC #18 with Diamond
No. 11307


Size / 160mm

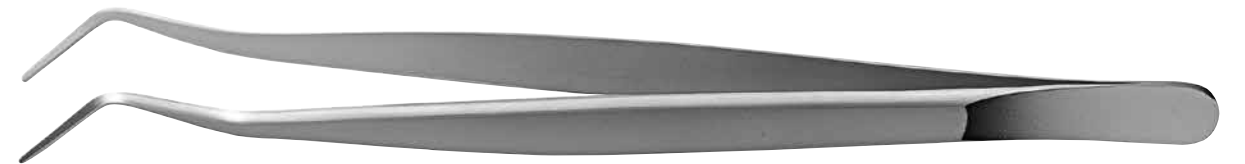

Premium
No. 11019

Size / 160mm

*Delicate tweezers.
*Careful finish prevents tip misalignment.

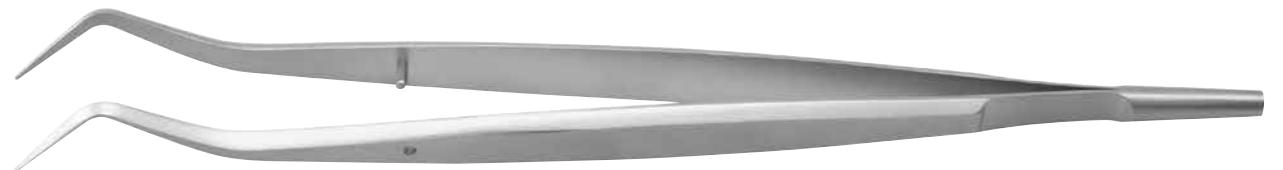

Koshi
No. 11020

Size / 165mm




Light Weight
No. 11024

Size / 165mm

-Clear 2-

*With thin glass, there is little blur of an image and reflected vividly.

*Suitable only for YDM mirror handle.

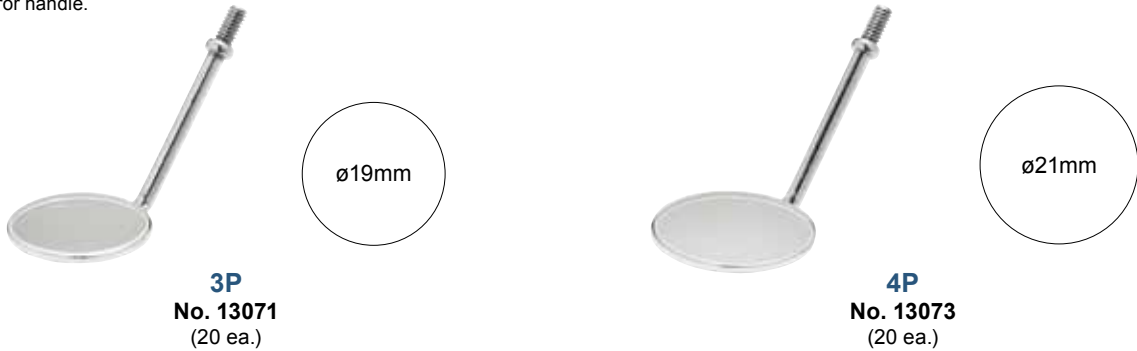

*Thickness / 1.9mm


*Reflection rate / About 93%

-Front Surface-

***Titanium Coating**

*Direct reflection prevents an image blurring.

*Suitable only for YDM mirror handle.

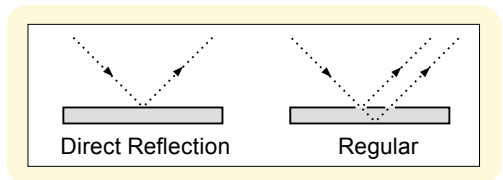


*Thickness / 1.9mm

*Reflection rate / About 98%

Metal Mirror (Simple Stem)

*The mirror made from stainless steel does not have a fear of being divided, and is suitable for the medical treatment of child dentistry, disabled person dentistry.

*Suitable only for YDM mirror handle.

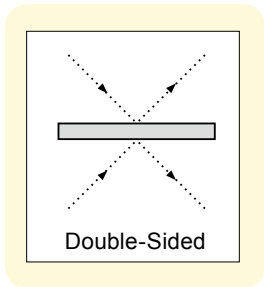

No. 13131
(6 ea.)




*Thickness / 1.1mm

*Reflection rate / About 65%

Mirror Handle

*Suitable only for YDM mouth mirror.

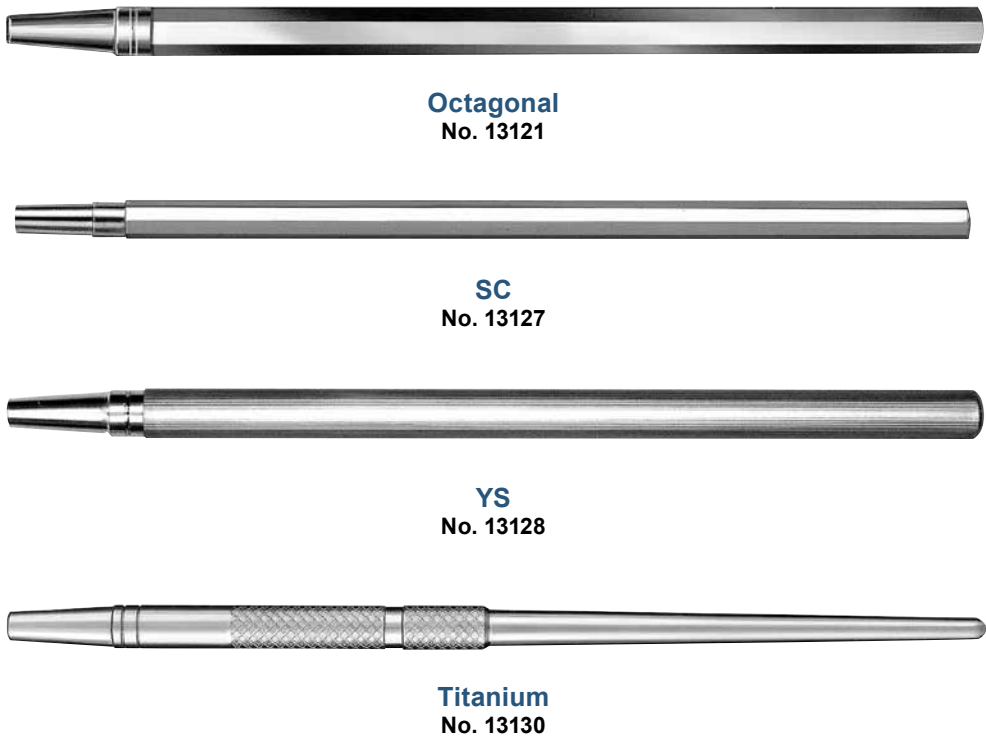

Mirror Handle with Millimeter Scale

*Suitable only for YDM mouth mirror.



*Material / POM

[How to measure endodontic file or reamer.]



①Fit stopper of reamer or file and pull in the direction of arrow.



②Measure the length of reamer or file with scale on the handle.



③Stopper of reamer or file is fixed.

*Measurement: Up to 30mm

Mouth Mirror
No. 28102
(50 ea.)

Explorer
No. 28101
(50 ea.)

Periodontal Probe
No. 28107
(50 ea.)

*Disposable type (Not sterilized)

*Material / PS : Heatproof temperature 70°C (Handle)

*Need disinfection by solution or EOG sterilization before use.

Temporary Filling Carrier

No. 13200

*A tip is heated and used.

*To prevent the burn of membrane, a black cover is made to slide to a tip.

Temporary Filling Remover

#10
No. 03010

Evacuation Tip

-Universal-
 *There is no trap for suction prevention attached at the evacuation tip, so it is easy to clean. (Rubber sleeve has a trap)

20°
No. 13611
(5 ea.)

90°
No. 13612
(5 ea.)

Sold Separately

Rubber Sleeve
No. 13613
(5 ea.)

*Size / 40mm

*Material / Si

***"Rubber Sleeve" has a trap inside.

Adapter
No. 13614

*Material / PTFE

*In case the inside diameter of evacuation tube is 16mm, this adaptor is needed for the use of "Evacuation Tip Universal".

Cleaning Brush L (ø12mm)
No. 13653
(12 ea.)

Size / 180mm

10

11



Stainless Dish



*Excellent when using small portion of bath such as topical anesthesia, local disinfection, etching solution, and temporary filling material (powder or liquid).

No. 13670
(5 ea.)



*Hole diameter / $\phi 16\text{mm}$
*Size / W56 x D31 x H8.5mm

Stainless Tray



M
No. 13413



Size / W245 x D146 x H8mm

S
No. 13410



Size / W245 x D73 x H8mm

