



Conservative Treatment

Excavator	38
Filling Instrument	40
Composite Syringe	44
Amalgam Product	44
Matrix Retainer / Separator	45
Contact Protector	45
Applicator	46
Rubber Dam Product	47
Endodontic Product	50
Fiber Cutter	58





Excavator



-Round Tip-




Str. #1 / Cvd. #1

Str. #2 / Cvd. #2



Str. #3 / Cvd. #3

Str. #4L / Cvd. #4L

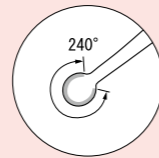*Material / Si (Handle)




*The edge along with the circle at a tip remove softened dentin and make cavity preparation efficiently.

Description		#1	#2	#3	#4L
Str.	 Cavity preparation for wall and floor.	01601	01602	01603	01604
Cvd.	 Used for undercut area, narrow rim of an opening and at the cavity angle.	01701	01702	01703	01704

-CDR-

*Excavator which pursued minimal intervention.

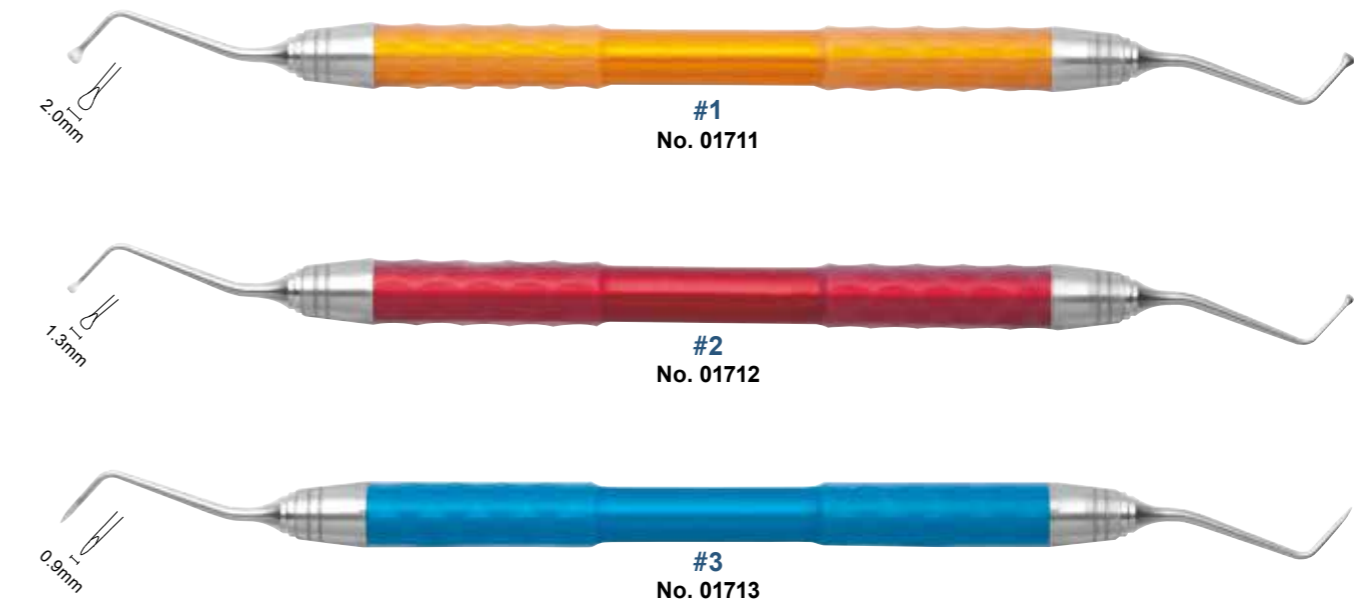


#1
No. 01711

#2
No. 01712

#3
No. 01713


*Material / Si (Handle)

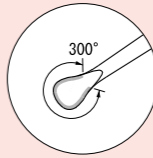

#1 / #2


*300° wide edge allows efficiently removal of cavity dentine.
*Slope working part design allows easy clean of particles.

Plastic Filling Instrument

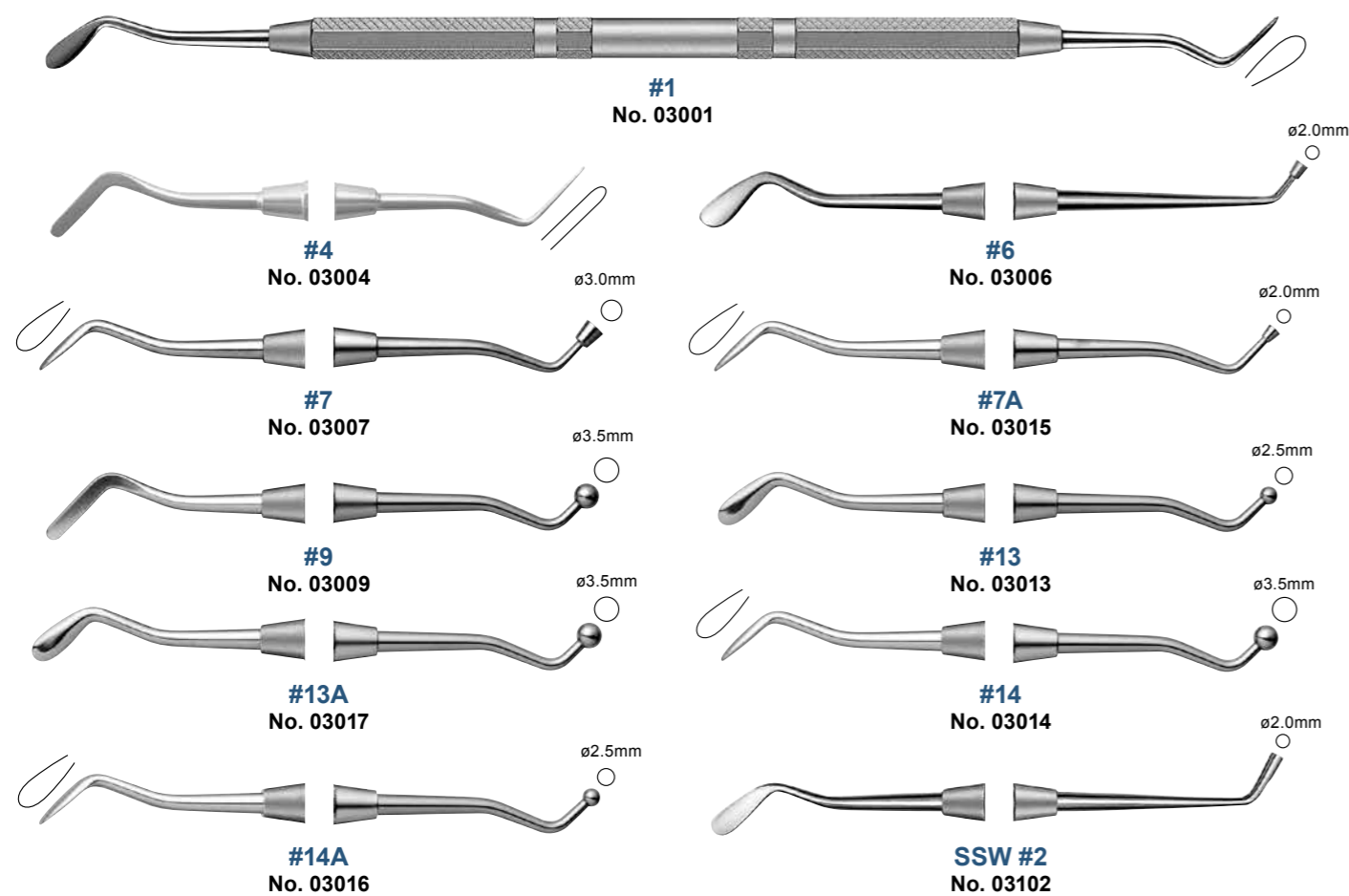

Proximal Filling Instrument

*For efficient filling to proximal, thickness at the tip is thin.

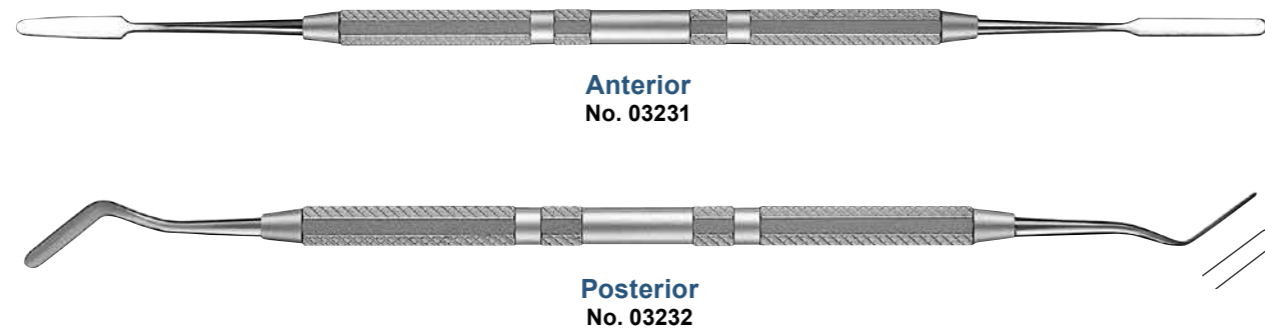

Cement Filling Instrument



No. 03241




Composite Filling Instrument

-DLC-

*Diamond Like Carbon Coating

*Remarkably stick-resistant for easy workability. Improved strength, and superior wear resistance.

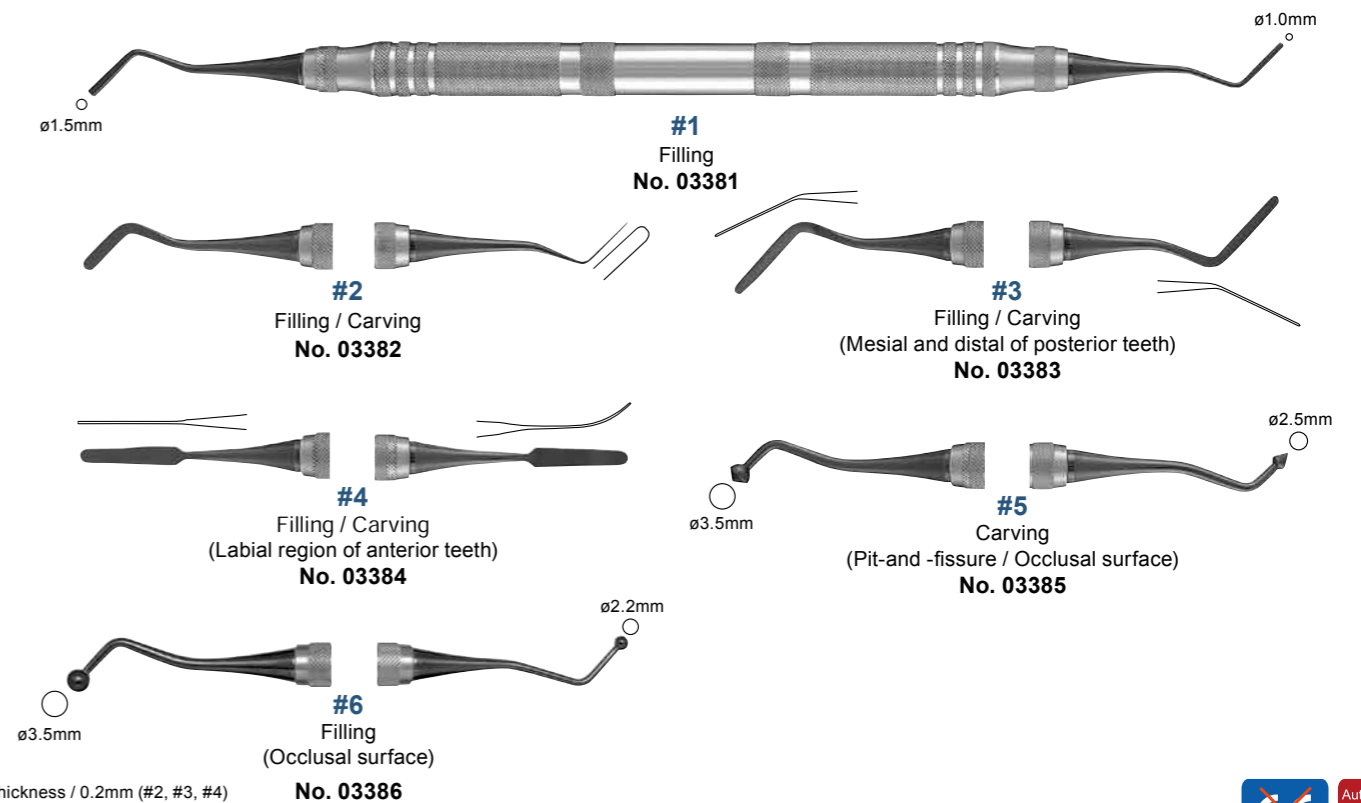


Composite Filling Instrument Tip (Replaceable Type)

*Remove all joint parts and replaceable tip before sterilization and disinfection.


*For more replaceable type products see P. 201.



-Handle-

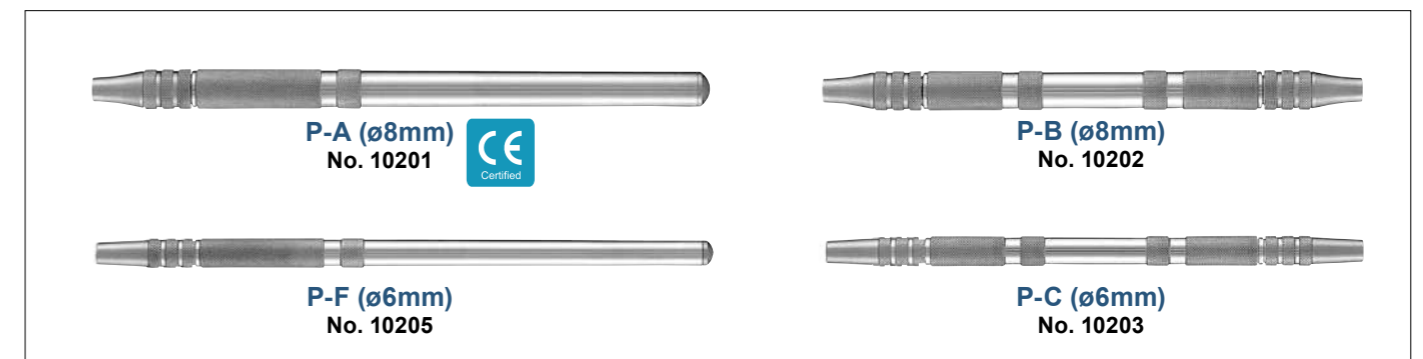

-Tip (3 ea.)-

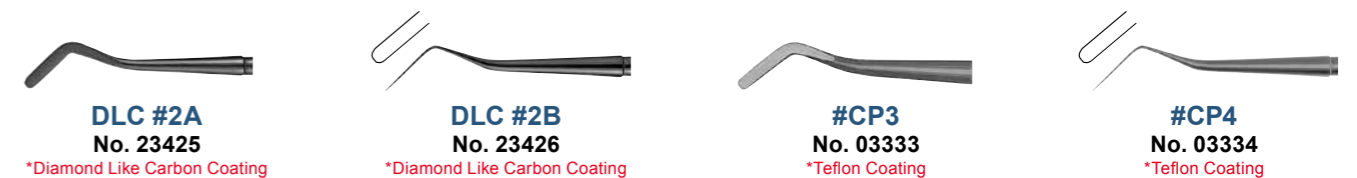

*Thickness / 0.2mm




Composite Filling Instrument






-TMDU 2 (Tokyo Medical Dental University)-

*Diamond Like Carbon Coating

*Remarkably stick-resistant for easy workability. Improved strength, and superior wear resistance.

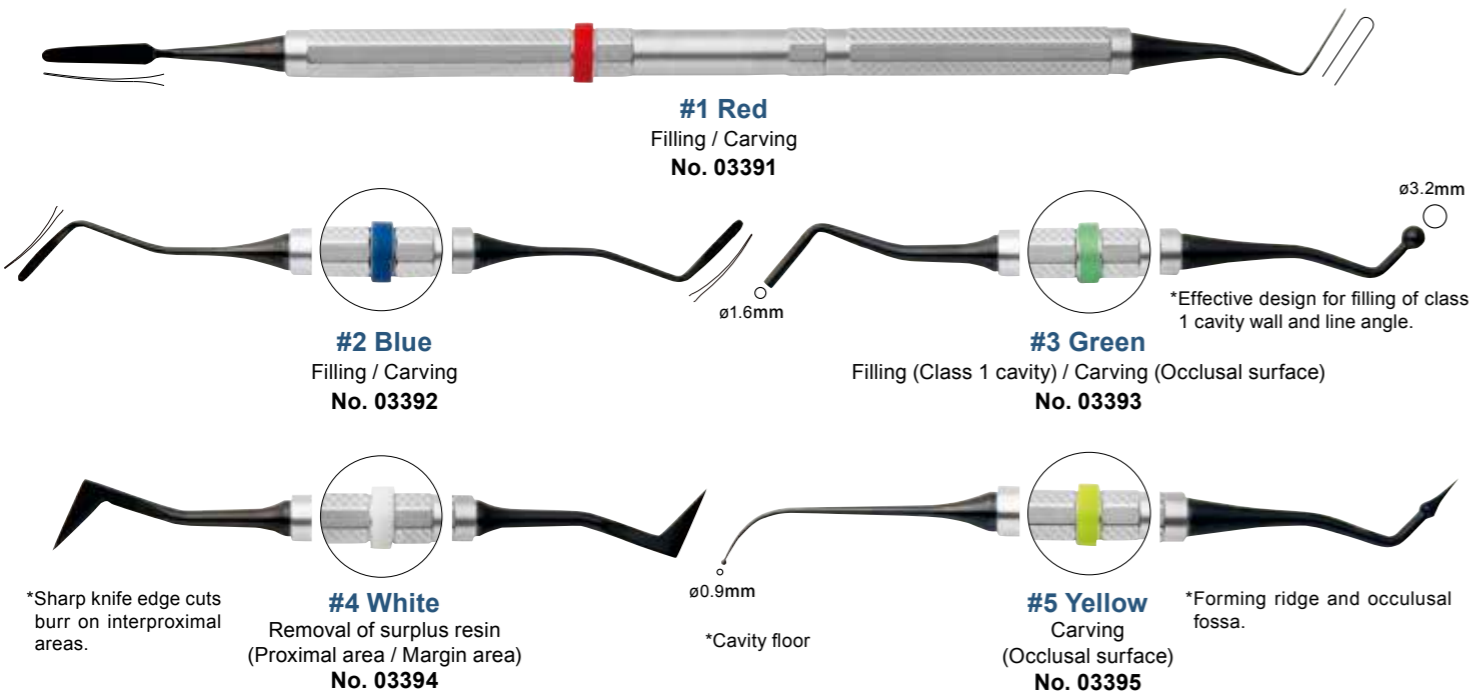





*Thickness / 0.2mm (#1)
0.1mm (#2)
*Material / Si (Color Ring)

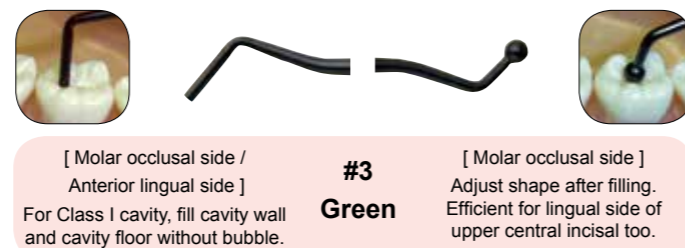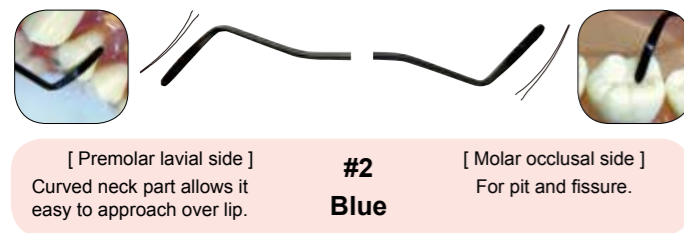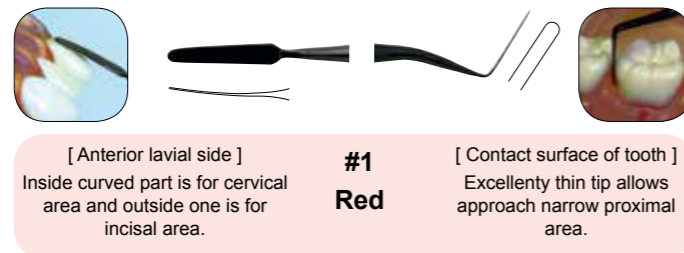

-TMDU 3 (Tokyo Medical Dental University)-


*Diamond Like Carbon Coating

*Remarkably stick-resistant for easy workability. Improved strength, and superior wear resistance.

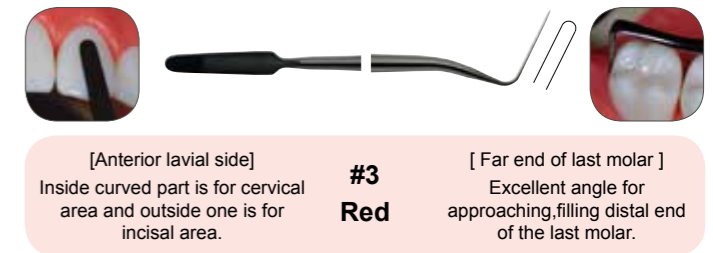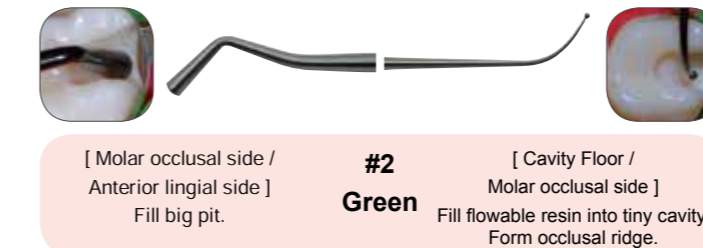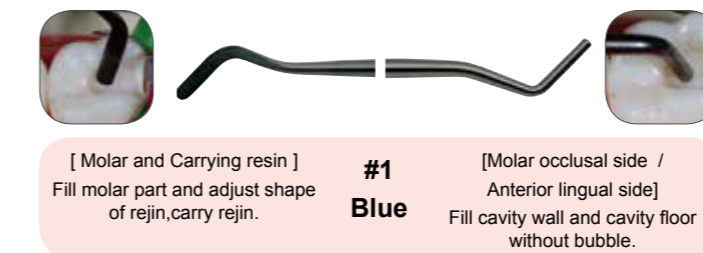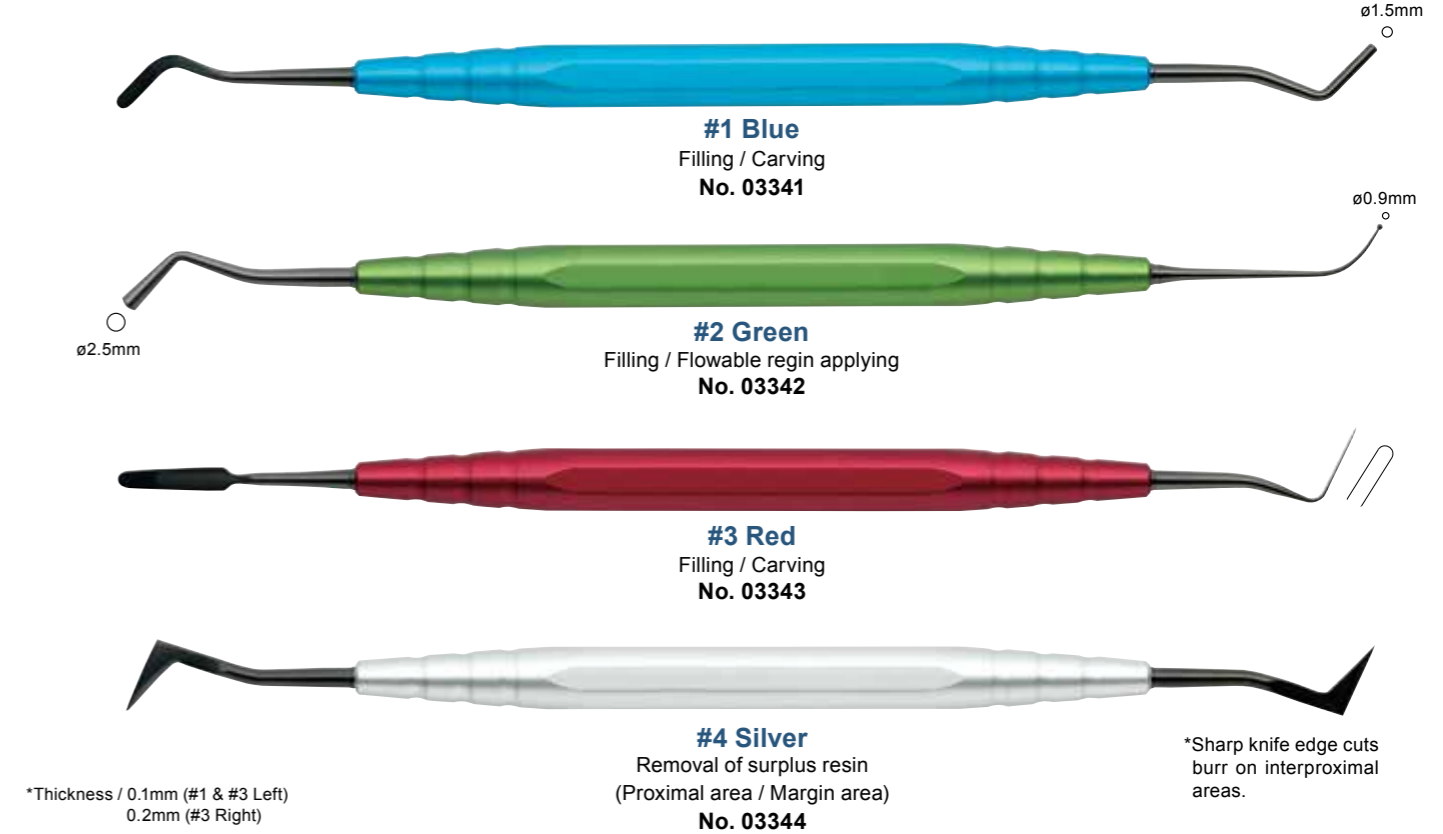

Contents

L : Main body L / Nozzle & Plug L (100 ea.)
S : Main body S / Nozzle & Plug S (100 ea.)

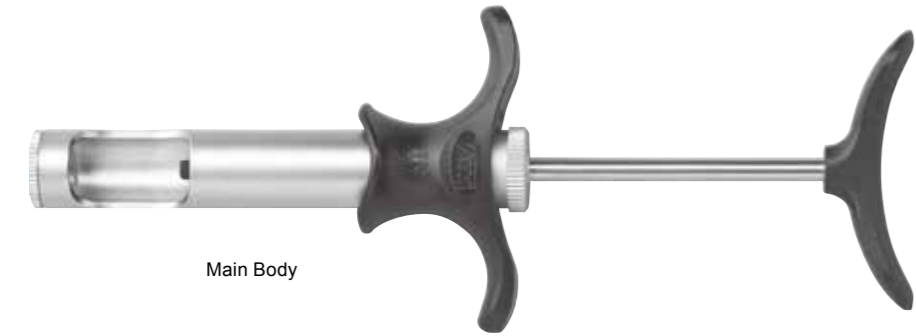

Main Body


Nozzle & Plug L
(Plug / Blue)



L

No. 13251


Nozzle & Plug S
(Plug / Green)


S

No. 13252

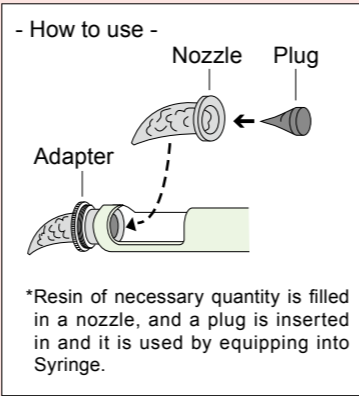

Individual Parts

No.	Description
13255	Main Body L
13256	Main Body S
13188	Nozzle & Plug L (100 ea.)
13190	Nozzle & Plug S (100 ea.)
13193	Nozzle & Plug L (500 ea.)
13194	Nozzle & Plug S (500 ea.)

*Material / PPS (Black part)

Amalgam Burnisher



S

No. 05002

Amalgam Carver



Round Handle

No. 07100


Matrix Retainer Set



Contents

Main body / Band (10 ea.)

No. 13170

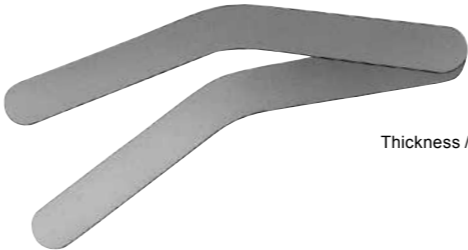

Thickness / 30μm

Individual Parts

No.	Description
13171	Main Body
13172	Band (10 ea.)

- How to use -


1) The both ends of a matrix retainer band are united and made a ring.

Crown part




Dental neck part

2) The matrix retainer band which became a ring is inserted in a matrix retainer and a fixed screw is tightened.

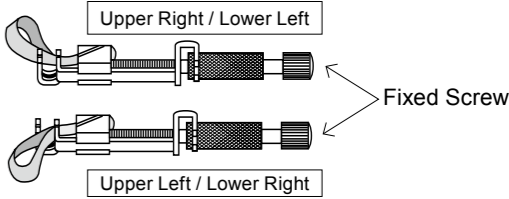

3) Apply to the target tooth and tightened adjustable screw.

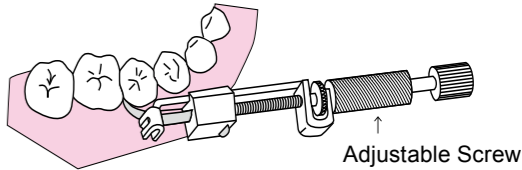

Separator



Ivory

Anterior

No. 13161



Elliot

Posterior

No. 13162



Contact Protector Set

Contents

Anterior (10 ea.) / Posterior (10 ea.)

*Stainless Steel


*Protect adjacent tooth during the cavity cutting.

No. 13175


Anterior



Posterior

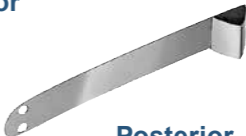

- How to use -

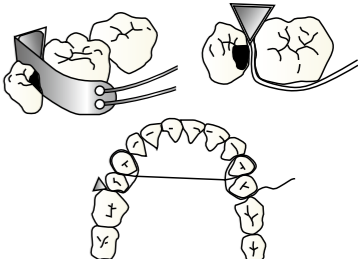

*It is used after letting froth pass in a hole.

*Thickness / 50μm

*Not sold separately.




Lining Applicator



-Round-



-Octagonal-

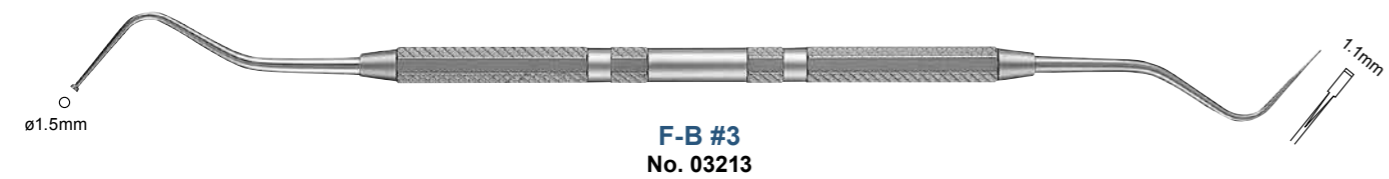

Ball-Pointed Applicator



No. 25015



Rubber Dam Clamp




***Stainless Steel**

*#12A / 13A → Serrated jaws and offset jaw height provides stability to partially erupted teeth.

*#26 / #26A → Used in pedodontic. In that case, #26 is for deciduous molars and #26A is for deciduous anterior teeth.



-Clamp shape-

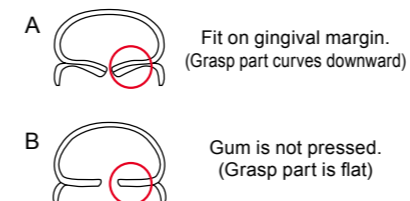


No.	Description		Shape
08526	#26	Upper molars	A
08501	#201	Upper molars	A
08556	#56	Lower molars	A
08502	#202	Lower molars	B
08506	#206	Premolars	A
08507	#207	Premolars	B
08527	#26A	Premolars	B

**Board
No. 08258**



Size / W132 x D87mm

- How to use -

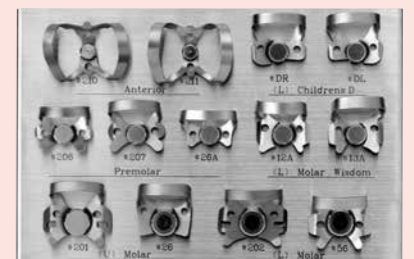

Reamer Guard

- *Can be arranged "Reamer Guard" by reamer or file type and size.
- *Cover will not allow the reamers and files to fall out.
- *An Engine reamer can be stored with "Bur Guard".

No.	Description
08431	Yellow
08432	Orange
08433	Pink
08434	Blue
08435	Green

-Holding length-

- *Storage capacity / 8 reamers and files
- *Size / W74 x D10 x H52 mm
- *Material / POM
- *Sterilizable with liquid sterilizing solutions.
- *Product does not include reamers and files.

Sold Separately

Case
No. 08439



- *Storage capacity / 18 ea.
- *Size / W167 x D126 x H62mm
- *Sterilizable with liquid sterilizing solutions.
- *Product does not include "Reamer Guard".



*The cover becomes a stand.



*Has a built in gauge to measure stop positions during root canal treatments.




*Able to connect each "Reamer Guards".

Reamer Stand

-Carry Set-

Contents Carry (3 ea.) / Case

No. 08411

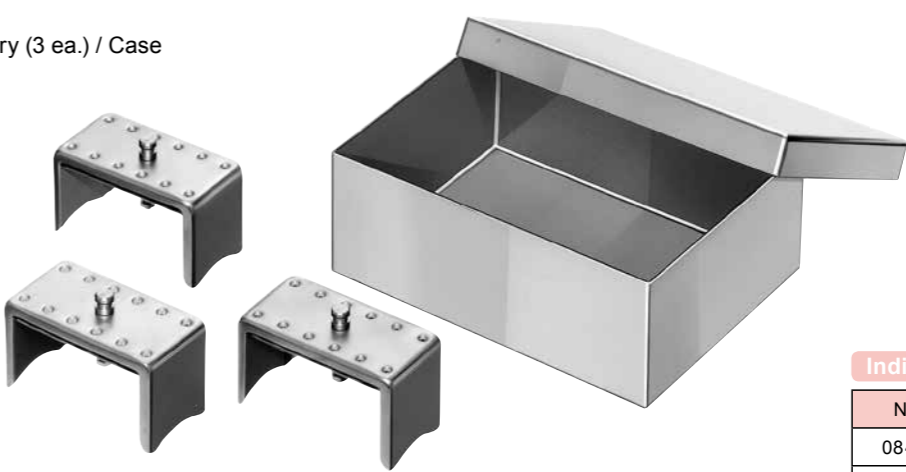

Individual Parts

No.	Description
08412	Carry (3 ea.)
24016	Case

- *Storage capacity / 12 reamers and files (Carry)
/ 6 ea. (Case)
- *Size / W60 x D30 x H33mm (Carry)
/ W120 x D90 x H40mm (Case)

-Safety-

- *It is ideal for managing reamers and files for every patients.
- *Can be autoclave or ultrasonic cleaning with reamers and files.
- *Lid does not allow reamers and files to fall out.

No. 08415



*Has a build in ruler.



- *Storage capacity / 11 reamers and files
- *Holding length / From 21 to 31mm
- *Size / W48 x D43 x H47mm
- *Material / PP (Green part)
- *Product does not include reamers and files.

Clean Pot

Contents Main Body / Paper (100 sheets)

No. 08423

- *During the treatment, can be cleaned reamers and files by paper.
- *Set 2 sheets to stabilize reamers for cleaning.



Individual Parts


Paper
No. 08422
(300 sheets)



- *Size / ø45 x H36mm
- *Product does not include reamers and files.



Contents Main Body / Joint (ø10mm & ø11mm Tapered) / Nozzle (5 ea.) / Cleaning Wire

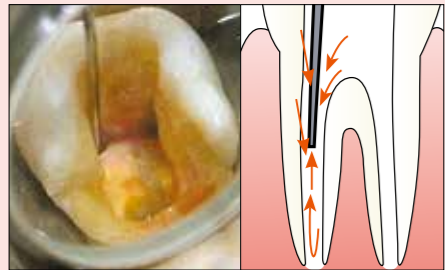

*The inside of a root canal is dried by an air current with suction in a root canal.
*Because of drying by an air current, it also becomes saving of a paper point.

*Tip outside diameter / ø0.6mm

Individual Parts



Sold Separately

Cleaning Brush S (ø4mm)
No. 13651
(12 ea.)



Size / 200mm

Adapter
No. 13614



*Material / PTFE

*In case the inside diameter of evacuation tube is 16mm, this adaptor is needed for the use of Evacuation Tip Universal.



Irrigating Syringe




3ml.
No. 13227

*Material / PFA (Cap / Main body)
PP (Trigger)

Explorer



for Root Canal Therapy #16
No. 01316

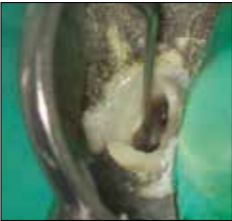

Broach Holder

No. 13150




File Holder Type-2

*Excellent for use in hard-to-reach posterior teeth.
*Applicable for root canal treatment of both mesial and distal.

No. 08064



Distal

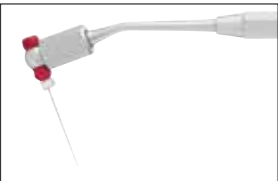

Mesial




Finger Ruler

No. 21155



Scale

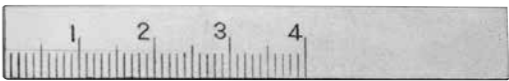


4cm
No. 21151



Endodontic Gauge

No. 21154

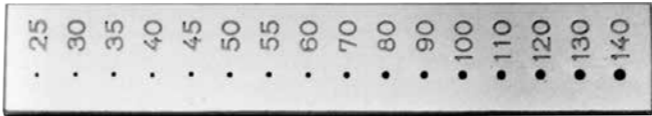

Front side





Back side



Root Forceps

*For removal of reamers and files deeper within root canal.
*The tip has a groove to securely hold residual reamer or file.

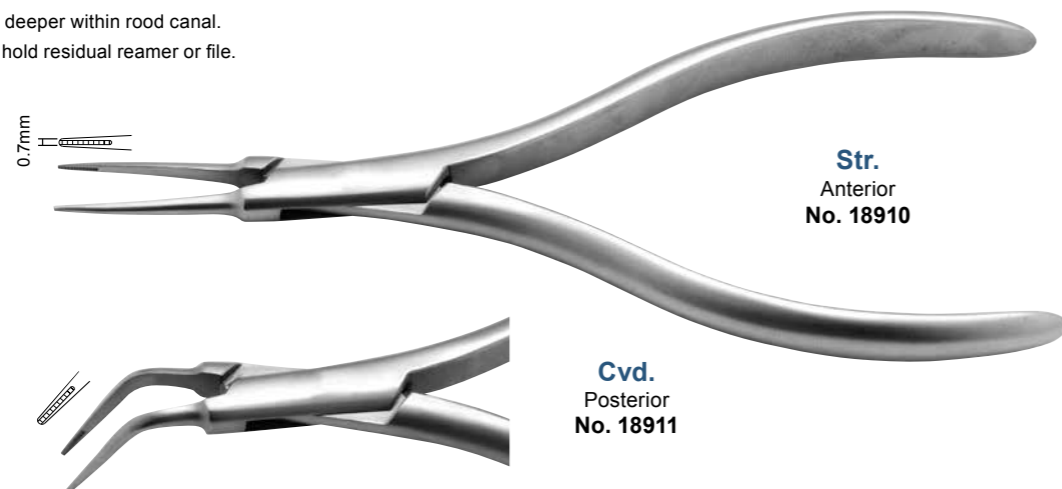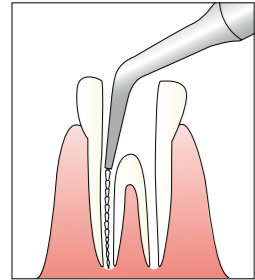


Str.
Anterior
No. 18910

Cvd.
Posterior
No. 18911

Stieglitz Root Forceps

*Assures firm grip of residual reamer or paper point in root canal.

No. 08063

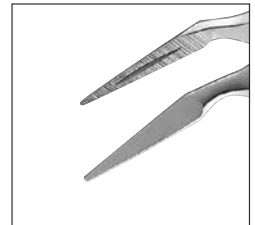

X-Ray Gauge


Contents Sleeve (2 ea.) / Gauge (2 ea.)

No. 21157

Sleeve



Gauge
(1mm by 1mm with hold line every 5mm)

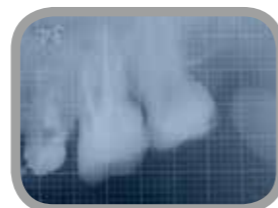

- How to use -



*Size / W38 x D48 x H8.5mm (Sleeve)
W33 x D45mm (Gauge)
*Material / Si (Sleeve)

*Suitable for root canal treatment, accurate examination before placing an implant, checking of calculus quantity and its area.

Gutta Percha Case

*Can be autoclave with points.

Mini
No. 24017



*Size / W91 x D42 x H13mm
*Material / Si (Green part)
*Product does not include main and accessory points.



Endodontic Tweezers




Size / 160mm



Size / 160mm





Size / 160mm

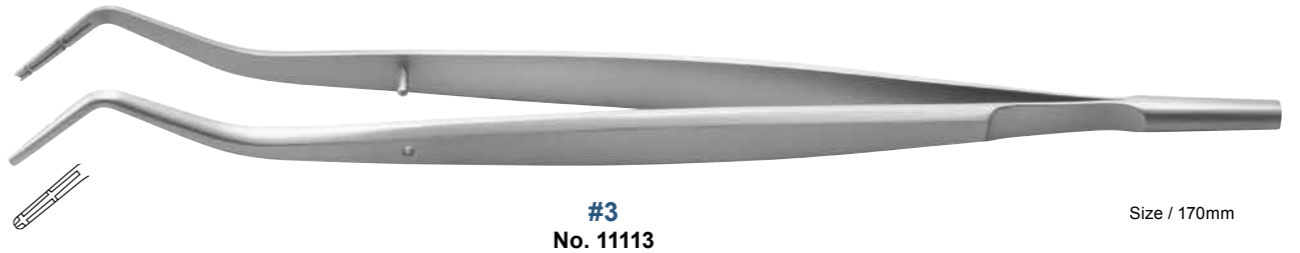

Size / 170mm



Size / 161mm

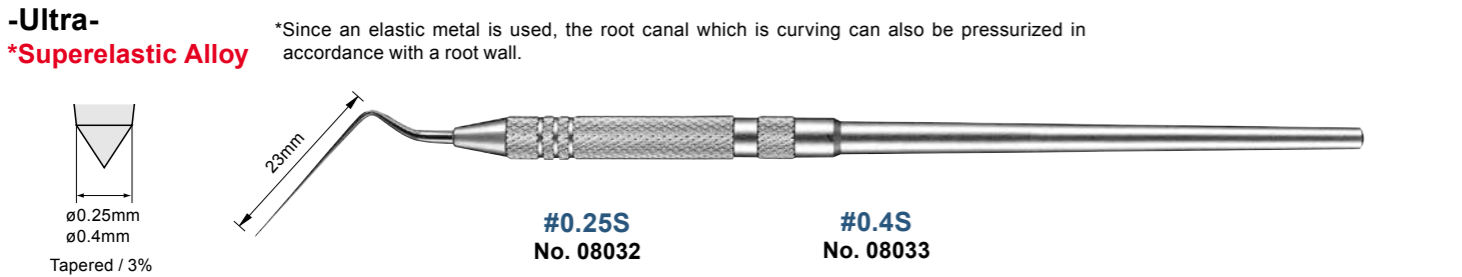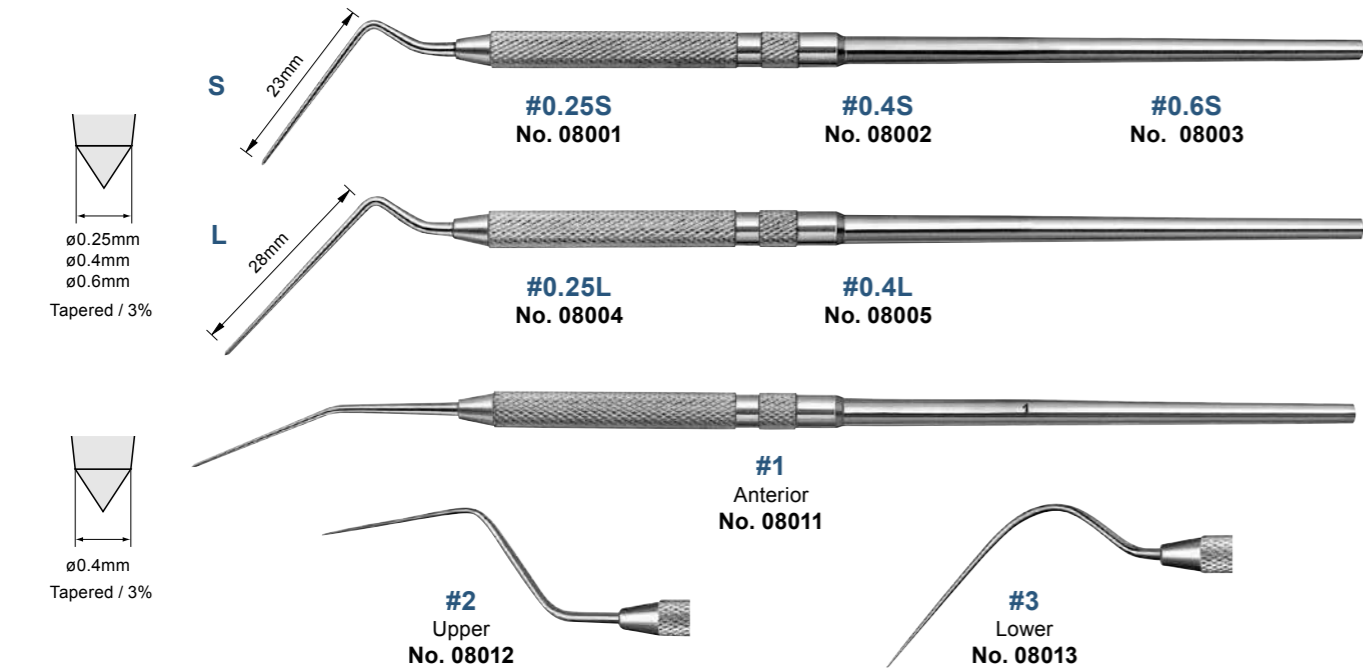



Fiber Cutter Type2

Contents Main body / Holder / Sleeve Guide (3 ea.)

- *Easy to cut the fiber post length for your need.
- *There is no fluff of a cutting plane.
- *Particulates are not generated.

No. 08071



- *Cut fiber can be entered in to the holder.
- *Easy to cut the fiber thanks to the gauge on the holder.
(Measurement from 5mm to 15mm)

*Size / 130mm
*Material / TPX (Holder)



Fiber Cutter Type2	Nipper	Steal Round Bur	Cutting Disc	Cutting Disc with Diamond
				

Problems from using Steal Round Bur, Cutting Disc, Diamond Disc

- 1) Fiber become fluffy.
- 2) Melt the surface because of the friction heat between fiber post and bur.
- 3) Remain the color of bur or disc on the section.